



Evolution[®] Esophageal

CONTROLLED-RELEASE STENT

Features

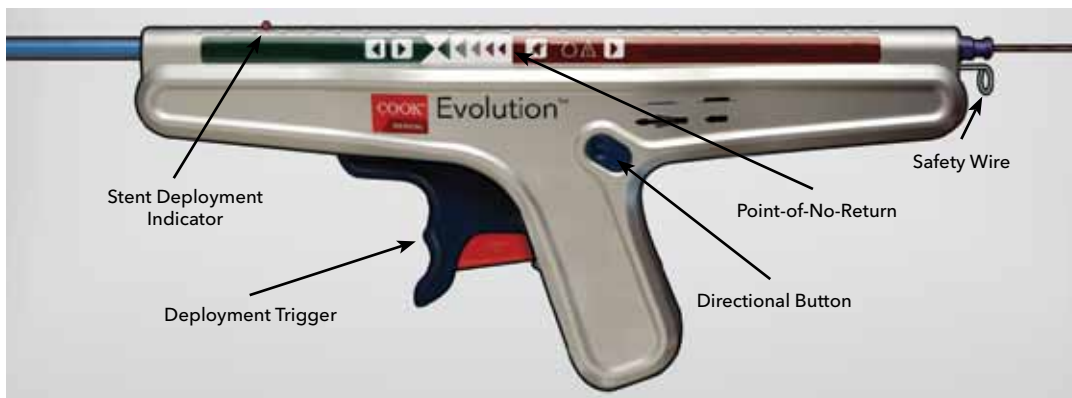


Step 1

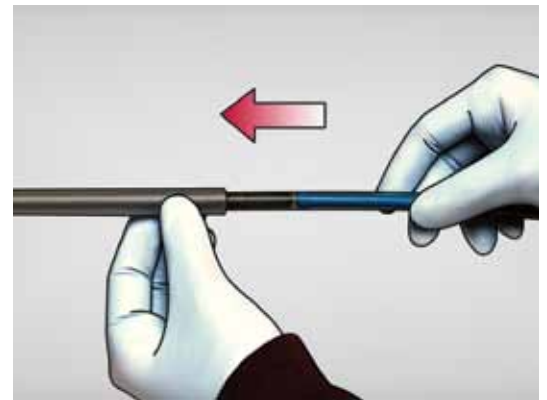

Remove protective tubing from stent.

Step 2

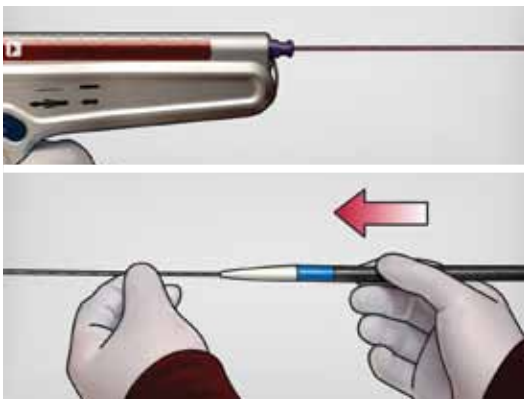

Using fluoroscopic guidance, introduce and advance introduction system in short increments over pre-positioned wire guide into esophagus.

Step 3



Fluoroscopically position radiopaque markers on inner catheter beyond extremities of stricture to be crossed.

Step 4



Confirm desired stent position fluoroscopically and deploy stent by removing red safety guard from handle. Continue deploying stent by squeezing trigger.


Step 5a



If repositioning is required during deployment, it is possible to recapture stent.

Note: It is not possible to recapture stent after passing point-of-no-return, indicated when the stent deployment indicator on top of introducer has passed point-of-no-return indicator on handle.


Step 5b

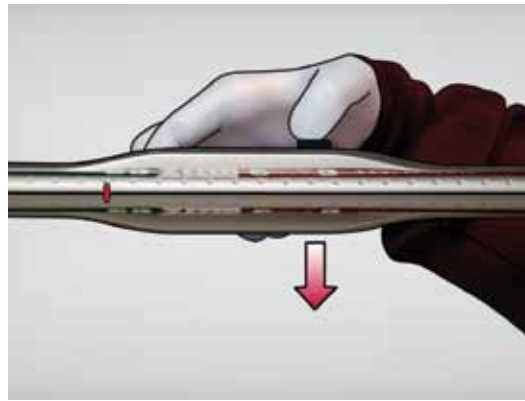

To recapture stent, push directional button on side of delivery system to opposite side.

Note: Hold thumb on button when squeezing trigger for first time to recapture. Continue squeezing trigger as required to recapture stent by desired amount.

Step 6

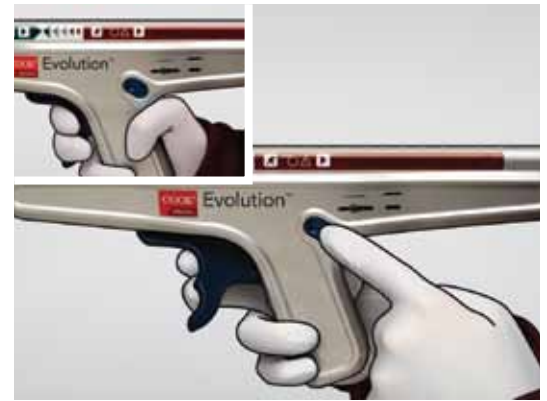


To resume deployment, push directional button to opposite side again and hold button for first stroke while squeezing trigger.

Step 7




When stent point-of-no-return has been passed, pull safety wire out of delivery handle near wire guide port. Continue deploying stent by squeezing trigger.

Step 8

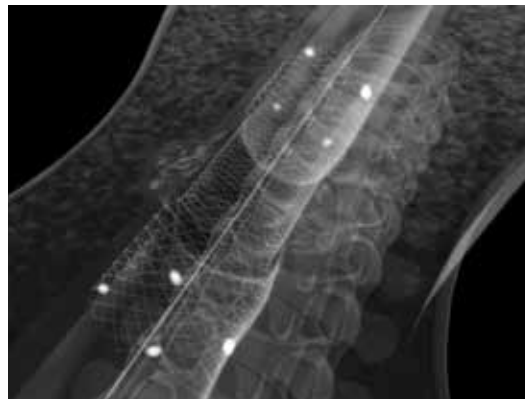

After deployment, fluoroscopically confirm stent expansion. Once expansion is confirmed, introduction system can be safely removed.

Step 9

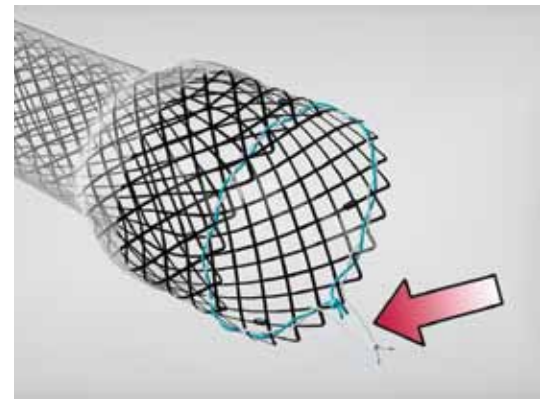

With the white lasso on proximal end of stent you can reposition stent, directly after placement by pulling with forceps on the white lasso and reposition stent to desired area.

Refer to current instructions for detailed system use.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do
Distributors: +353 61239240, ssc.distributors@cookmedical.com
Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001633, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Norway: +47 23162968, no.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com

www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
E-mail: orders@cookmedical.com

Australia:
Phone: +61 738411188, 1800777222, Fax: +61 738411288, 1800077283
E-mail: cau.custserv@cookmedical.com

