
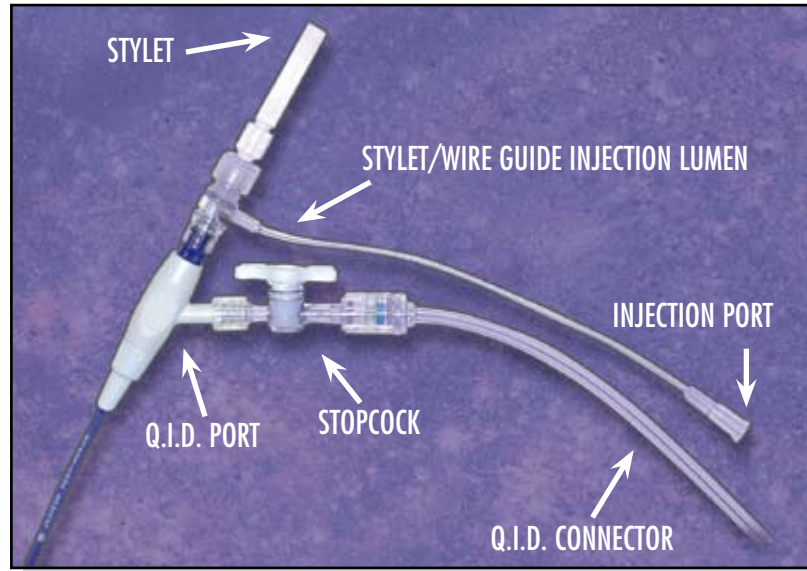


QUICK REFERENCE GUIDE

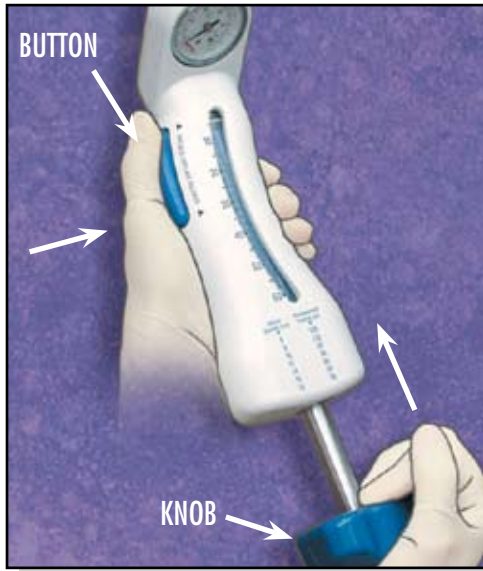

STEP 1

Pressurize/Depressurize
Ensure the pressurize/depressurize knob is rotated completely counter-clockwise and the plunger is pushed completely forward. To move the plunger, depress and hold the inflate/deflate button on the side of the Q. I. D. Release the button to lock the plunger in place.

STEP 2

Fill Inflation Device
Submerge the connecting tube in a minimum of 60 cc water or saline. Fill the inflation device with a minimum of 40 cc and a maximum of 50 cc of fluid by depressing the button while simultaneously pulling back on the knob. Release the button and lock the plunger in place.



STEP 3

Purge Air
To completely purge the inflation device of air, hold the Q.I.D. upright and depress the button while simultaneously grasping and slowly pushing forward on the knob. Release the button to lock the plunger in position.



STEP 4

Attach Dilator
Attach the Luer Lock fitting of the Q.I.D. connection tube to the side port of the y-connecting cap of the balloon dilator. If use of the stopcock is desired, attach the appropriate fittings to the balloon dilator and the Q.I.D. connecting tube. The arm of the stopcock indicates which port is closed.




STEP 5

Maintain Vacuum
Depress the button and pull back on the knob until vacuum is indicated on the pressure gauge. While maintaining a vacuum with the knob, release the button. Release the knob. The set-up of the Q. I. D is now complete.

QUICK REFERENCE GUIDE

STEP 6

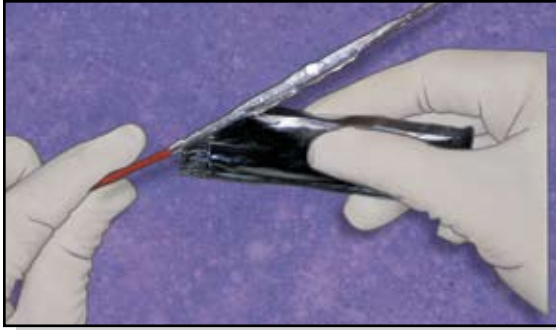

Introduction and Positioning

Prior to introduction into the endoscope, lubricate balloon. **Do not pre-inflate balloon.** Prior to insertion of balloon dilator, negative pressure is mandatory to maintain balloon deflation. Advance balloon to the desired position. **NOTE:** Radiopaque metal bands located at either end of the balloon may be used to fluoroscopically confirm balloon position. **NOTE:** To use a wire guide, remove the stylet and flush the lumen. Insert wire guide.





STEP 7

Balloon Inflation

For rapid inflation, depress button and push knob forward until significant resistance is encountered. Maintain pressure on the knob while releasing button.

Rotate the knob clockwise and increase to the specified balloon pressure. See Catheter Tag for balloon pressure information.

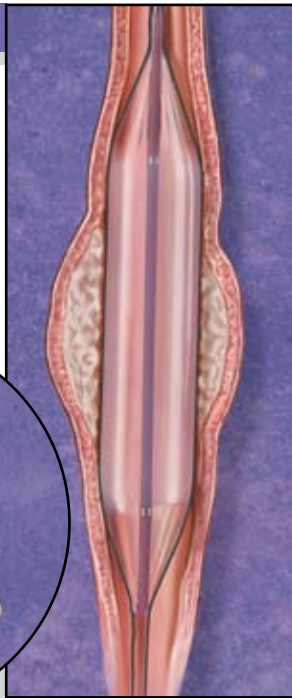
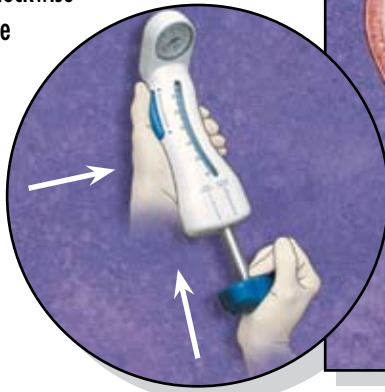


STEP 8

Balloon Deflation

Decrease balloon pressure by rotating the knob counterclockwise until pressure decreases below 20 psi. Depress button and pull back on the knob, creating negative pressure as seen on the pressure gauge. Release button to lock in negative pressure. The balloon can now be removed.

NOTE: If the balloon does not appear fully deflated, you may need to purge air from the inflation device. Disconnect the Q.I.D. from the balloon. Hold the Q.I.D. upright and depress the button while simultaneously grasping and slowly pushing forward on the knob. Release the button to lock the plunger in position. If the optional stopcock is attached, open stopcock to atmosphere and purge inflation device of air by closing the port to the balloon catheter. Depress the button and push forward on the knob.

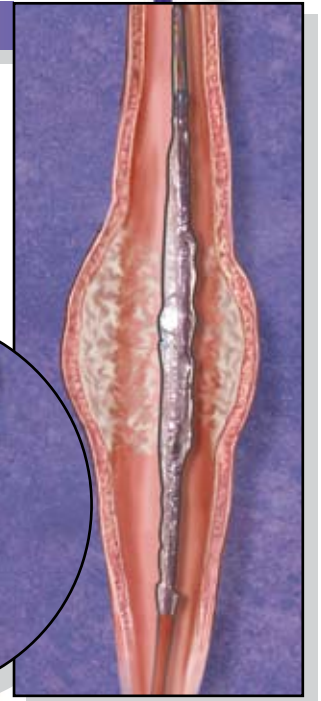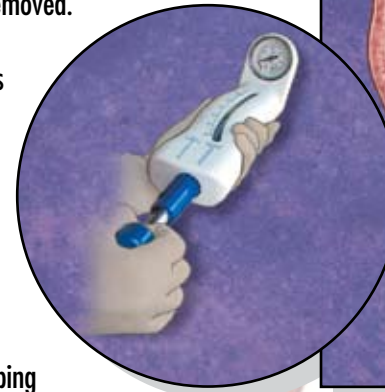

CATHETER TAGS



COOK[®]

Endoscopy

4900 Bethania Station Rd • Winston-Salem, NC 27105
 336-744-0157 • Fax: 800-743-1147 • 336-744-1147
USA Customer Service: 800-457-4500

www.cookendoscopy.com

Refer to current instructions for detailed system use.

© 2006 Cook Endoscopy.