


COOK
MEDICAL

Biodesign™
SURGISIS® HERNIA GRAFT

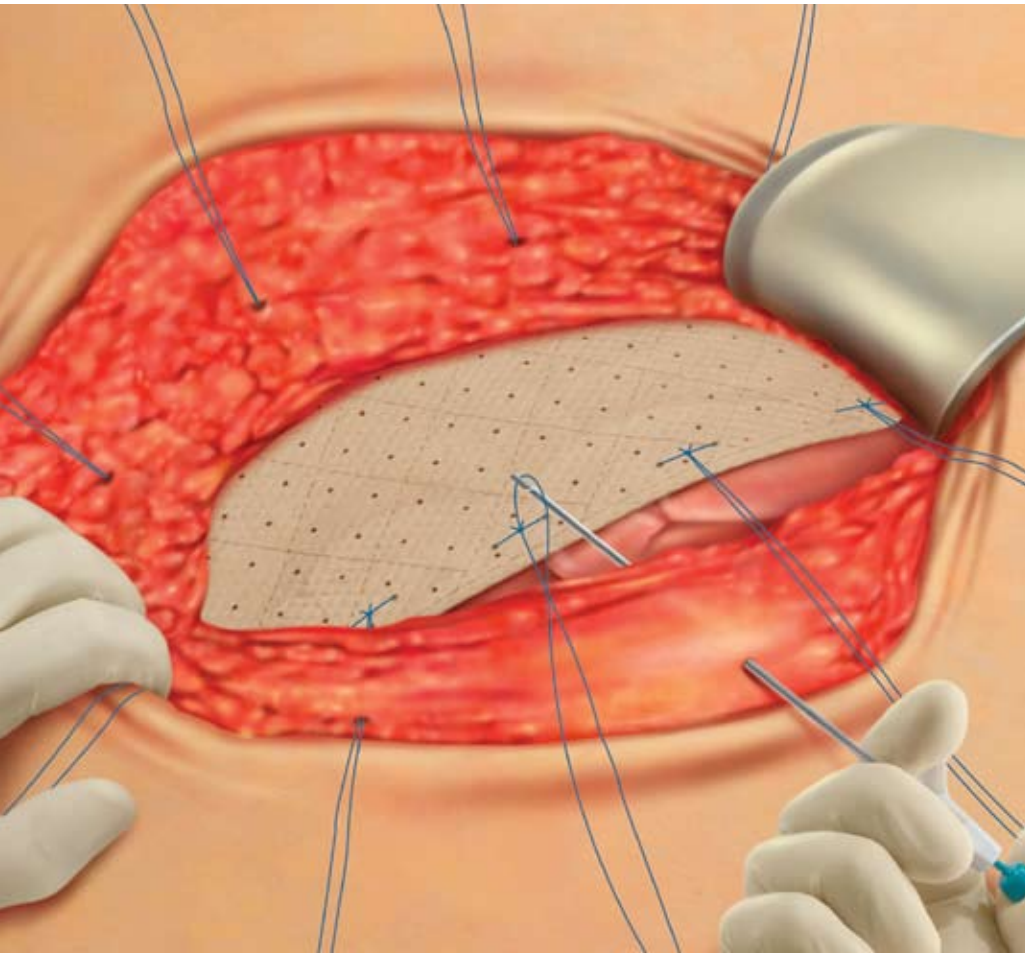

Restore long-term strength.


The Biodesign Hernia Graft helps you start fast and finish strong when repairing hernias. Biodesign is a scaffold that lends support while helping damaged or infected tissue restore itself quickly.

Biodesign signals the body to promote rapid and complete remodeling, is resistant to infection and encapsulation, and offers reduced early complications. Over time, it becomes strong, fully vascularized tissue that creates a lasting repair.

As a tension-free repair method, Biodesign provides a more stable solution and quicker patient recovery than primary repair methods. It rehydrates quickly, can be used in either laparoscopic or open surgery, and is available in large sizes.

Cook Medical—Advancing patient care through biologic design.

Global Product Number	Order Number	Size cm
Hernia Graft		
G36032	C-SLH-8H-13X15	13 x 15
G46600	C-SLH-8H-13X22	13 x 22
G36033	C-SLH-8H-20X20	20 x 20
G48216	C-SLH-8H-20X30	20 x 30

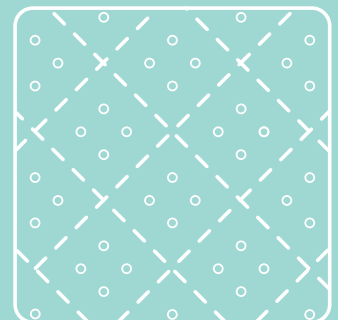

FEATURES & BENEFITS

- Supports tissue as it remodels; restores natural form and function.
- Fast-hydrating material can be ready to suture in less than 1 minute.
- Provides long-term strength with minimal stretching and tension-free repair.
- Open-matrix design is easy to handle and allows rapid remodeling.

TIPS

- For best results, underlay **Biodesign 5 cm beyond edges of hernia**.
- **Use only antiseptics and antimicrobials that have** proven compatible with wound healing. Problems have been reported with gentamicin, povidone iodine, bacitracin, polymyxin B and vancomycin.
- **Slow-resorbing sutures** are recommended in conjunction with an **interrupted suture technique**.
- **Space sutures 2 to 3 cm apart**.
- **Reapproximating tissue** is recommended for adequate graft-to-tissue contact.

DATA

1. Badylak S, Kokini K, Tullius B, et al. Strength over time of a resorbable bioscaffold for body wall repair in a dog model. *J Surg Res.* 2001;99(2):282-287.
2. Badylak S, Kokini K, Tullius B, et al. Morphologic study of small intestinal submucosa as a body wall repair device. *J Surg Res.* 2002;103(2):190-202.
3. Badylak SF, Park K, Peppas N, et al. Marrow-derived cells populate scaffolds composed of xenogeneic extracellular matrix. *Exp Hematol.* 2001;29(11):1310-1318.
4. Badylak SF. Small intestinal submucosa (SIS): a biomaterial conducive to smart tissue remodeling. In: Bell E, ed. *Tissue Engineering: Current Perspectives*. Cambridge, MA: Burkhauser Publishers; 1993:179-189.
5. Edelman DS. Laparoscopic herniorrhaphy with porcine small intestinal submucosa: a preliminary study. *JSLs.* 2002;6(3):203-205.
6. Franklin ME Jr, Gonzalez JJ Jr, Glass JL. Use of porcine small intestinal submucosa as a prosthetic device for laparoscopic repair of hernias in contaminated fields: 2-year follow-up. *Hernia.* 2004;8(3):186-189.
7. Gupta A, Zahriya K, Mullens PL, et al. Ventral herniorrhaphy: experience with two different biosynthetic mesh materials, Surgisis and Alloderm. *Hernia.* 2006;10(5):419-425.
8. Helton WS, Fisichella PM, Berger R, et al. Short-term outcomes with small intestinal submucosa for ventral abdominal hernia. *Arch Surg.* 2005;140(6):549-560.
9. Hodde JP, Hiles MC. Bioactive FGF-2 in sterilized extracellular matrix. *Wounds.* 2001;13(5):195-201.
10. Kropp BP, Sawyer BD, Shannon HE, et al. Characterization of small intestinal submucosa regenerated canine detrusor: assessment of reinnervation, in vitro compliance and contractility. *J Urol.* 1996;156(2, pt 2):599-607.
11. Rehydration comparative studies. Data on file, Cook Biotech, Inc.
12. Angiogenesis and rapid incorporation studies. Data on file, Cook Biotech, Inc.
13. Franklin ME Jr, Trevino JM, Portillo G, et al. The use of porcine small intestinal submucosa as a prosthetic material for laparoscopic hernia repair in infected and potentially contaminated fields: long-term follow-up. *Surg Endosc.* 2008;22(9):1941-1946.

PROCEDURAL AREAS

HERNIA REPAIR

FISTULA REPAIR
PLASTIC & RECONSTRUCTIVE/ENT
STAPLE LINE REINFORCEMENT
PEYRONIE'S REPAIR
CONTINENCE RESTORATION
DURAL REPAIR
PELVIC FLOOR REPAIR

VENTRAL/INCISIONAL HERNIA

HERNIA GRAFT
INGUINAL HERNIA
HIATAL HERNIA
PARASTOMAL HERNIA
ABDOMINAL LOCK

EVOLUTION OF TISSUE REPAIR



Biodesign is a breakthrough technology that provides signals and support for the body to restore itself. It incorporates the best attributes of a biologic graft—complete remodeling and resistance to infection—yet it is easy to use, widely available and moderately priced.



www.cookmedical.com

COOK MEDICAL INCORPORATED

P.O. Box 4195, Bloomington, IN 47402-4195 U.S.A.
Phone: 812.339.2235, Toll Free: 800.457.4500, Toll Free Fax: 800.554.8335

COOK (CANADA) INC.

111 Sandiford Drive, Stouffville, Ontario, L4A 7X5 CANADA
Phone: 905.640.7110, Toll Free: 800.668.0300

WILLIAM A. COOK AUSTRALIA PTY. LTD.

95 Brandl Street, Brisbane Technology Park, Eight Mile Plains
Brisbane, QLD 4113 AUSTRALIA, Phone: +61 7 3841 1188

COOK IRELAND LTD.

O'Halloran Road, National Technology Park, Limerick, IRELAND
Phone: 353 613 34440

© COOK 2009 SUR-BM-EXL-EN-200904

AORTIC
INTERVENTION

CARDIOLOGY

CRITICAL
CARE

ENDOSCOPY

INTERVENTIONAL
RADIOLOGY

PERIPHERAL
INTERVENTION

SURGERY

UROLOGY

WOMEN'S
HEALTH