

This changes everything.



StoneBreaker™

SALIVARY PNEUMATIC LITHOTRIPTER

Ground-breaking

The StoneBreaker™ Salivary Pneumatic Lithotripter eliminates the need for open surgery to remove salivary stones, expanding endoscopic interventional solutions in the salivary ducts.

Powerful


StoneBreaker uses pneumatic force to quickly and directly fragment salivary stones in the parotid and/or submandibular ducts.

Cost-effective

StoneBreaker is a reusable, portable and nonelectric interventional tool.



www.cookmedical.com

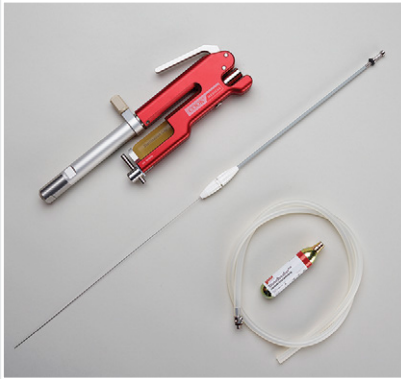

StoneBreaker™

SALIVARY PNEUMATIC LITHOTRIPTER

The StoneBreaker Salivary Pneumatic Lithotripter is a portable, compact, nonelectric intracorporeal lithotripter intended to fragment stones in the parotid and/or submandibular ducts. StoneBreaker is powered by a CO₂ cartridge and requires no electrical or pneumatic connections. The mechanical shock that the device generates is transferred as kinetic energy through the salivary probe to the tip, which is placed directly against the stone that is to be fragmented.

StoneBreaker Salivary Pneumatic Lithotripter and components

Reusable Products

Order Number	Reference Part Number	Description
G34852	SSM-1	StoneBreaker Salivary Pneumatic Lithotripter
G34926	SSMA-SC	Sterilization Cap
G34927	SSMA-PC	Salivary Probe Cap
G34928	SSMA-EL	Exhaust Line
G34929	SSMA-EC	Exhaust Cap

Single-Use Products

Order Number	Reference Part Number	Description	Length mm	Working Length mm	Working Diameter mm
G34851	SBPA-056-500-255	StoneBreaker Salivary Pneumatic Lithotripter Probe	500	252	0.56
G34930	SSMC-10	CO ₂ Cartridge	-	-	-

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001633, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: orders@cookmedical.com

Australia:
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
 E-mail: cau.custserv@cookmedical.com

