


Simplifying Hemostasis.



ONE CLIP INFINITE POSSIBILITIES

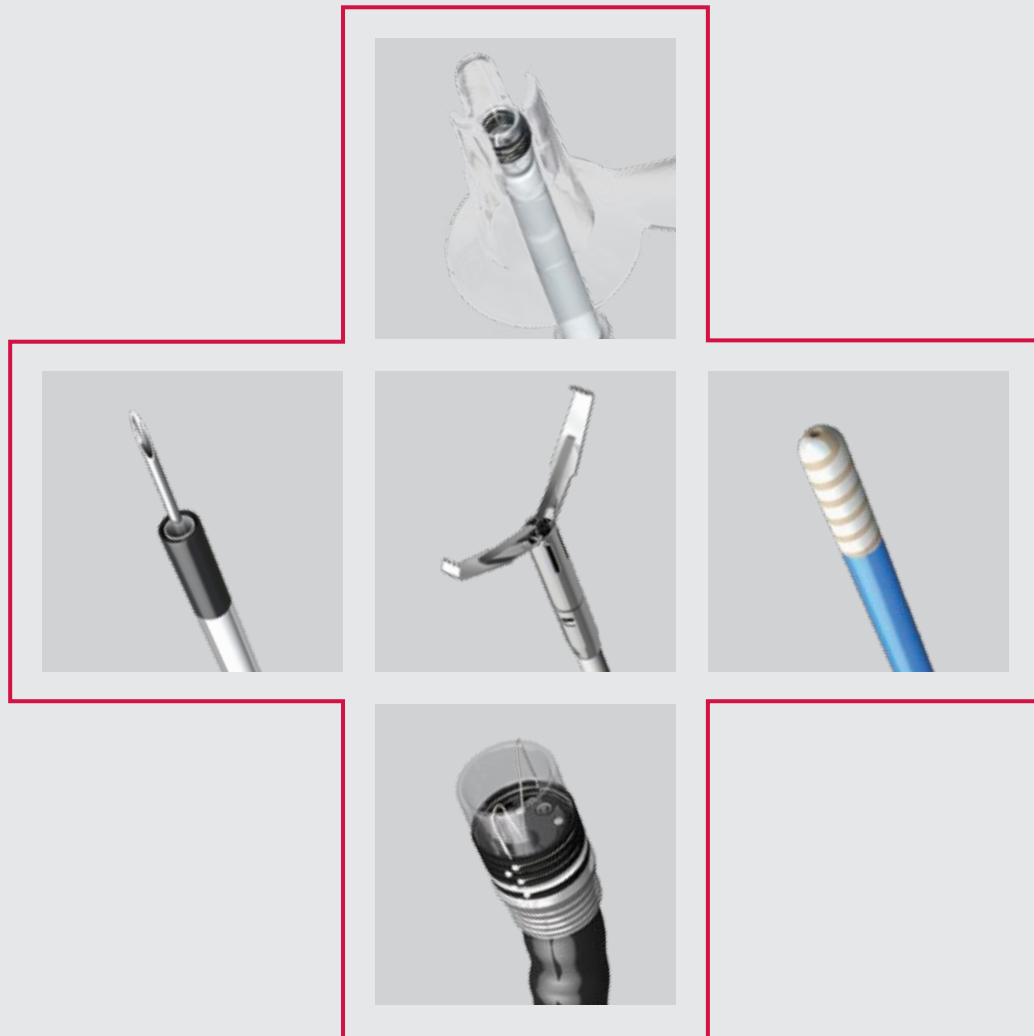

WIDEST SPAN



STRENGTH & SECURITY




360° ROTATION



MECHANICAL



Instinct™ Endoscopic Hemoclip

- Widest span on the market and fully adjustable so you can securely grasp as much, or as little, tissue as your procedure requires
- 360-degree, bidirectional rotation and open/close ability for precise clip placement
- Distinctive anchoring tips improve your tissue gripping capability while the robust, nitinol-reinforced jaws give you added stability
- Instinctive handle and clip design for simple, one-step clip detachment
- MR Conditional per ASTM F2503 



MECHANICAL



4,6,10 Shooter®
Saeed Multi-band Ligator



- 19 consecutive years as the market leader in ligation
- Opti-Vu® Barrel with recessed bands and friction fit adaptor designed to optimize your endoscopic view while maintaining suction and secure fit on endoscope
- Universal-fit barrels streamline order information and simplify inventory control

MECHANICAL



ShortShot®
Saeed Hemorrhoidal Multi-band Ligator with TriView™

- Preloaded for multiple hemorrhoidal ligations in a single procedure
- TriView anoscope with positioning slots to help optimize anatomical positioning of hemorrhoids
- Only a single insertion of an anoscope is required, which potentially decreases patient discomfort and procedure time

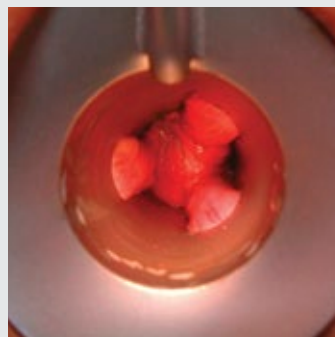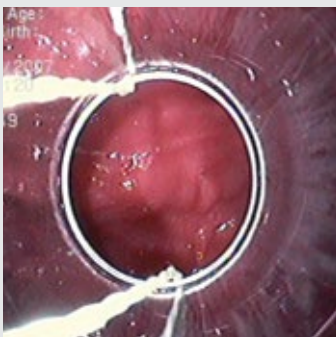

INJECTION



AcuJect®
Variable Injection Needle

- Smooth needle advancement and retraction
- Clear outer catheter with white inner catheter allows visualization of needle advancement
- Inner and outer catheters lock together to secure the advanced needle during injection
- Adjustment capability for desired needle length extension


THERMAL



Quicksilver®
Bipolar Probe

- Hydromer grafted tip and spiral electrodes designed to minimize sticking for effective coagulation
- Central lumen allows for irrigation
- 350 cm catheter with single or dual plug options for optimal scope and generator compatibility
- Completely disposable, eliminating cleaning and reprocessing costs



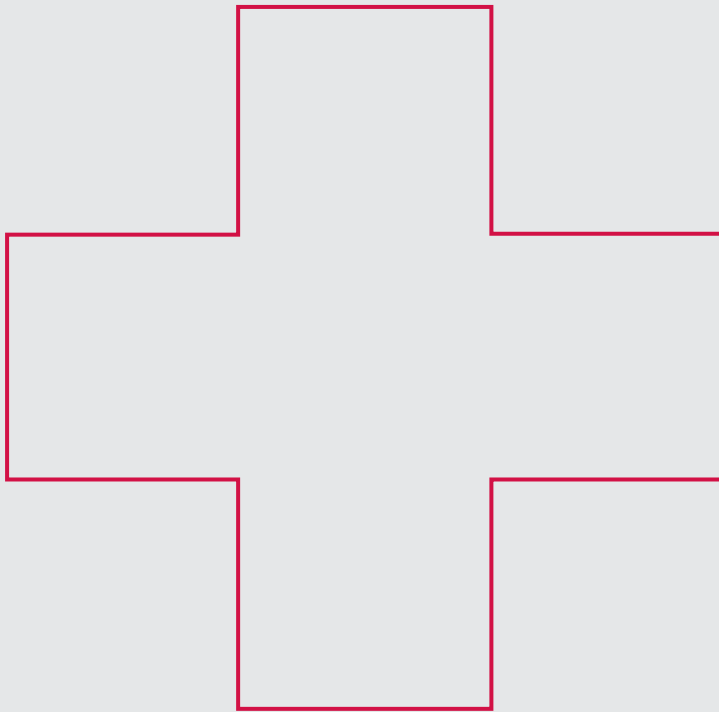

Simple solutions for complex problems.

When facing any GI bleed, the more simple and effective treatment options you have, the better your patient outcomes can be. That's why, for the past 30 years, Cook has given you a wide variety of easy-to-use, efficient hemostasis options.

From the intuitive control of the Instinct Endoscopic Hemoclip and the market-leading Saeed Six Shooter Multiband Ligator to a selection of time-tested needles and probes—see how our simple solutions can solve your complex problems.

Instinct™

Endoscopic Hemoclip

Order Number	Reference Part Number	Catheter Size Fr	Catheter Length cm	Box Quantity
G18343	INSC-7-230-S	7	230	10
G18344	INSC-7-230-20-S*	7	230	20

*Check for availability.

4,6,10 Shooter®

Saeed Multi-band Ligator

Order Number	Reference Part Number	Number of Bands	Endoscope Outer Diameter mm	Trigger Cord Length cm
G31915	MBL-U-4	4	8.6-11.3	122
G31917	MBL-U-6	6	8.6-11.3	122
G31919	MBL-U-10	10	8.6-11.3	122
G31916	MBL-U-4-F	4	8.6-11.3	142
G31918	MBL-U-6-F	6	8.6-11.3	142
G31920	MBL-U-10-F	10	8.6-11.3	142
G25119	MBL-6-XL	6	11-14	122
G25221	MBL-6-XL-C	6	11-14	190

ShortShot®

Saeed Hemorrhoidal Multi-band Ligator with Triview™

Order Number	Reference Part Number	Number of Bands
G34063	HMBL-4-TRI	4

Supplied with a TriView Anoscope.

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

AcuJect®

Variable Injection Needle

Order Number	Reference Part Number	Catheter Size Fr	Needle gage	Catheter Length cm
G22525	VIN-23	7	23	220
G22526	VIN-25	7	25	220

Quicksilver®

Bipolar Probe

Order Number	Reference Part Number	Catheter Size Fr	Catheter Length cm	Generator Plug Type
G24925	BCP-7A	7	350	dual
G24926	BCP-7B	7	350	single
G24927	BCP-10A	10	350	dual
G24928	BCP-10B	10	350	single

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001633, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: orders@cookmedical.com

Australia:
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
 E-mail: cau.custserv@cookmedical.com

