

From cytology to histology: Advance your EBUS diagnostic procedures.



EchoTip ProCore[®]

ENDOBONCHIAL HD ULTRASOUND BIOPSY NEEDLE

At Cook, we were among the first to recognize the clinical impact of EBUS procedures. That's why we focused our expertise as market leaders in ultrasound needles to create a full range of EchoTip ProCore EBUS needles to meet your endobronchial ultrasound needs. Now, for the first time, you can potentially achieve cytology and histology in your EBUS procedures, even in the challenging pulmonary hilar region. Our goal is simple: Give you the histological information to help you deliver the best patient care.




cookmedical.com




Olympus



Pentax

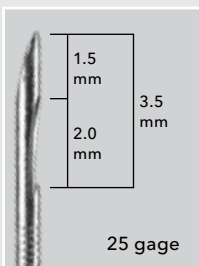

25 gage

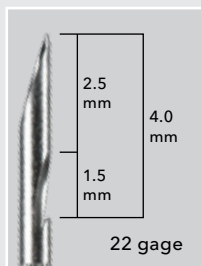

22 gage

EchoTip ProCore®

ENDOBONCHIAL HD ULTRASOUND BIOPSY NEEDLE

Olympus - Used with an ultrasound endoscope for fine needle biopsy, (FNB), of submucosal and extramural lesions within or adjacent to the tracheobronchial tree or gastrointestinal tract.

Pentax - Used with an ultrasound endoscope for fine needle biopsy, (FNB), of submucosal and extramural lesions within or adjacent to the tracheobronchial tree.

Supplied sterile and is disposable - intended for single use only.

- Featuring the core trap designed for receiving tissue in the delicate pulmonary area
- High Definition dimpling pattern for needle visibility
- Naturally contoured handle offers needle control & stability
- Adjustable sheath length
- 5 cm needle extension presents a large diagnostic reach
- On scopes with metal luer hubs, the removable sheath adjuster allows direct connection to the scope
- On scopes with metal non-luer hubs, the sturdy scope adaptor design allows:
 - secure connection and procedural stability
 - suction and injection capability due to the airtight valve

| Order Number | Reference Part Number | Needle Size Gage | Sheath Size Fr | Adjustable Needle Extension cm | Stylet Tip | Minimum Accessory Channel mm |
|--------------|-----------------------|------------------|----------------|--------------------------------|-------------------|------------------------------|
| G34281 | ECHO-HD-22-EBUS-O-C* | 22 | 4.1 | 0-5 | recessed ball tip | 2 |
| G34279 | ECHO-HD-22-EBUS-P-C** | 22 | 4.1 | 0-5 | recessed ball tip | 2 |
| G34282 | ECHO-HD-25-EBUS-O-C* | 25 | 4.1 | 0-5 | beveled | 2 |
| G34280 | ECHO-HD-25-EBUS-P-C** | 25 | 4.1 | 0-5 | beveled | 2 |

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

*Needle attaches to scope adapter on Olympus Endobronchial Ultrasound scope with non-luer metal hub.

**Needle attaches to luer metal hub on Pentax Endobronchial Ultrasound scope.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001633, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: orders@cookmedical.com

Australia:
 Phone: +61 738411188, 1800777222, Fax: +61 738411288, 1800077283
 E-mail: cau.custserv@cookmedical.com

