



Precision and control.

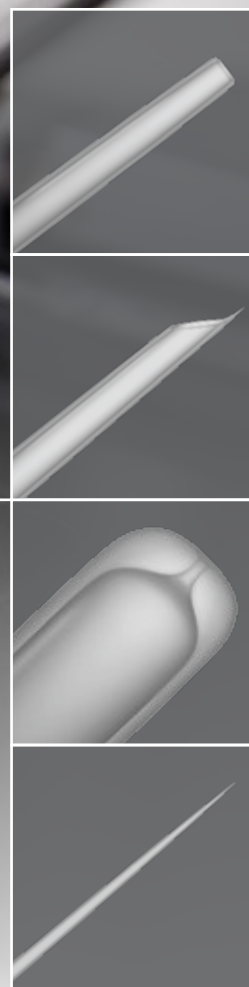

MICROMANIPULATION PIPETTES

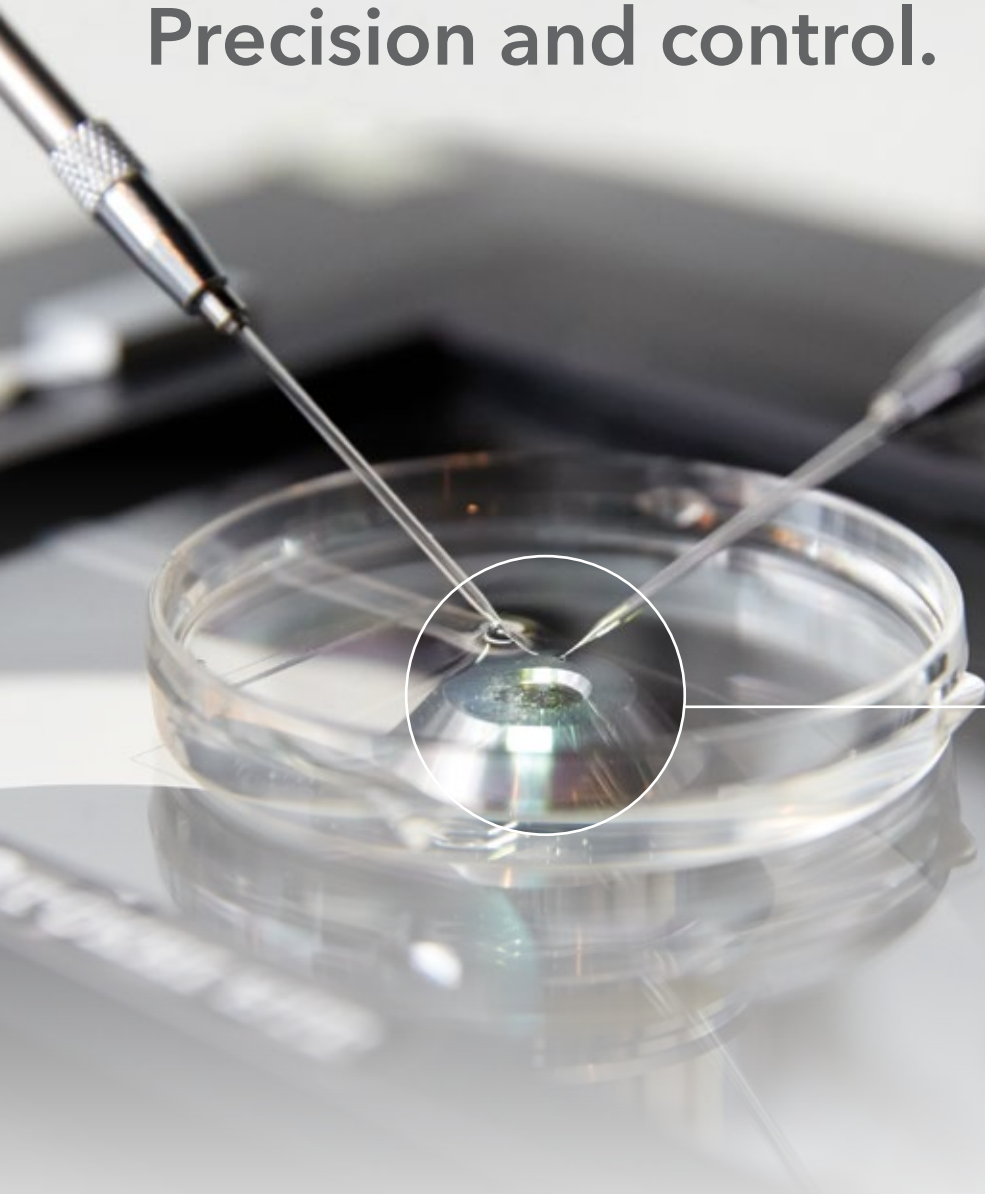

Precise

Cook Medical micromanipulation tools are precision instruments manufactured to exacting quality standards.

Specialized

Pipettes are available in different configurations to meet the needs of intracytoplasmic sperm injection (ICSI) and preimplantation genetic diagnosis (PGD) procedures.

Tested

All pipettes are inspected under extreme magnification, sterilized by gamma irradiation, and lot tested for one-cell MEA survivability and LAL toxicity.



Protected

Pipettes are packaged individually in containers that are designed to deliver undamaged product. These containers are easy to use and are made of recyclable material.

Cook® Microinjection Pipette

The microinjection pipette is used for the intracytoplasmic single sperm injection of oocytes. These pipettes have well-defined bevels to puncture the zona pellucida, and the tip-to-bend length is the same for all sizes.

Sold in boxes of 10.

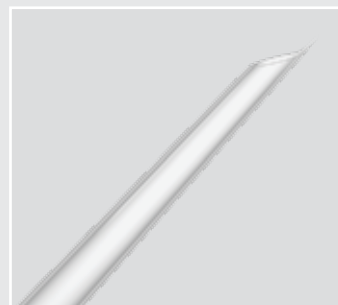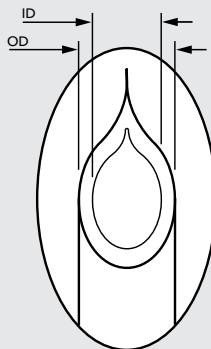


| Order Number | Reference Part Number | Inner Diameter μm | Outer Diameter μm | Distal Tip Angle degrees | Bevel Angle degrees | Tip-to-Bend Length μm |
|--------------|-----------------------|------------------------------|------------------------------|--------------------------|---------------------|----------------------------------|
| G32817 | K-MPIP-3130 | 4.7 | 6 | 30 | 30 | 700 |
| G32818 | K-MPIP-3135 | 4.7 | 6 | 35 | 30 | 700 |
| G32814 | K-MPIP-1030 | 5.0 | 7 | 30 | 30 | 700 |
| G18090 | K-MPIP-1035 | 5.0 | 7 | 35 | 30 | 700 |
| G32821 | K-MPIP-3330 | 5.5 | 7 | 30 | 30 | 700 |
| G26684 | K-MPIP-3335 | 5.5 | 7 | 35 | 30 | 700 |

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Cook® Holding Pipette

The holding pipette is used to hold an oocyte, embryo, or blastocyst in position with the application of vacuum during intracytoplasmic single sperm injection or assisted hatching/zona drilling. Featuring a smooth, flame-polished surface, the holding pipette provides stability and holding power.

Sold in boxes of 10.


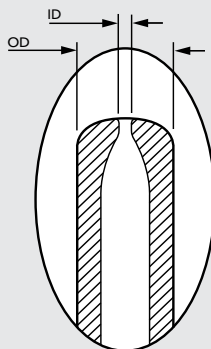


| Order Number | Reference Part Number | Inner Diameter μm | Outer Diameter μm | Distal Tip Angle degrees | Tip-to-Bend Length μm |
|--------------|-----------------------|------------------------------|------------------------------|--------------------------|----------------------------------|
| G32801 | K-HPIP-1030 | 17 | 80 | 30 | 700 |
| G32802 | K-HPIP-1035 | 17 | 80 | 35 | 700 |
| G32806 | K-HPIP-2130 | 23 | 110 | 30 | 700 |
| G32807 | K-HPIP-2135 | 23 | 110 | 35 | 700 |
| G32810 | K-HPIP-3330 | 34 | 120 | 30 | 700 |
| G32812 | K-HPIP-3335 | 34 | 120 | 35 | 700 |

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Cook® Blastomere Aspiration Pipette

Blastomere aspiration pipettes are used to aspirate a blastomere or trophectoderm to diagnose genetic disorders prior to embryo selection. The pipette tip is cut straight and flame-polished smooth.

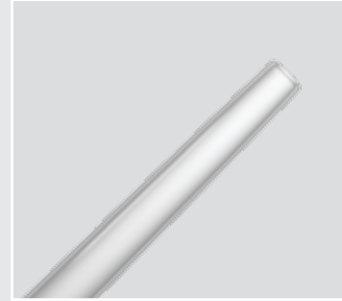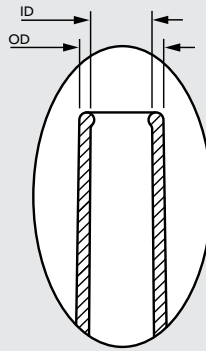



Sold in boxes of 10.

| Order Number | Reference Part Number | Inner Diameter μm | Outer Diameter μm | Distal Tip Angle degrees | Tip-to-Bend Length μm |
|--------------|-----------------------|------------------------------|------------------------------|--------------------------|----------------------------------|
| G32795 | K-EBPH-3035 | 30 | 42 | 35 | 2,500 |
| G32799 | K-EBPH-3535 | 35 | 49 | 35 | 2,500 |

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Cook® Adjustable Handle

The Cook Adjustable Handle is used with Flexipet pipettes for handling gametes and embryos. The aspiration volume can be preset at 0.25 μL to 3.0 μL . The handle accepts all sizes of Flexipet pipettes and provides the fine control that is needed to manipulate oocytes, embryos, and blastomeres. The accessory set consists of replacement parts for the handle, including one small o-ring, one large o-ring, a spacer, and a collet.

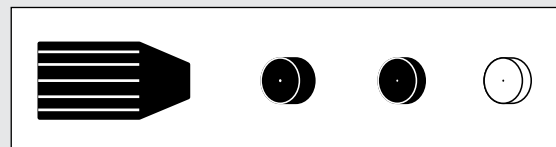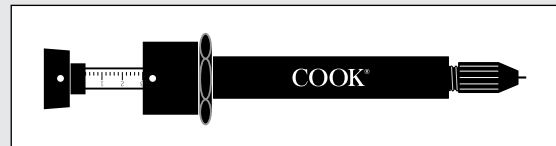

| Order Number | Reference Part Number |
|----------------------|-----------------------|
| Handle | |
| G18674 | K-MPH-1000 |
| Accessory Set | |
| G18904 | K-FRO-1000 |

Cook® Flexipet® Pipette

Flexipet Denuding Pipettes are used for blastomere or polar body manipulation, oocyte and embryo manipulation, or denuding.

Flexipet Manipulation Pipettes are used for blastomere or polar body manipulation, oocyte and embryo manipulation, blastocyst handling, and manipulation of the oocyte-cumulus complex.

- Made of flexible polycarbonate
- Designed to flex rather than scratch surfaces or break
- Can be used with the Cook Adjustable Handle



| Order Number | Reference Part Number | Inner Diameter µm | Application |
|---------------------------------------|-----------------------|-------------------|--|
| Flexipet Denuding Pipettes | | | |
| 5 Vials of 10 (50 Total) | | | |
| G46021 | K-FPIP-1120-10BS-5 | 120 | oocyte and embryo manipulation |
| G26711 | K-FPIP-1130-10BS-5 | 130 | oocyte and embryo manipulation |
| G26712 | K-FPIP-1140-10BS-5 | 140 | oocyte and embryo manipulation |
| G26713 | K-FPIP-1170-10BS-5 | 170 | oocyte and embryo manipulation |
| 10 Vials of 1 (10 Total) | | | |
| G46019 | K-FPIP-1140-10 | 140 | oocyte and embryo manipulation |
| G46018 | K-FPIP-1170-10 | 170 | oocyte and embryo manipulation |
| Flexipet Manipulation Pipettes | | | |
| 1 Vial of 10 (10 Total) | | | |
| G26181 | K-FPIP-1080-10BS | 80 | blastomere or polar body manipulation |
| G26718 | K-FPIP-1100-10BS | 100 | oocyte and embryo manipulation |
| G18878 | K-FPIP-1300-10BS | 300 | blastocyst handling |
| 5 Vials of 10 (50 Total) | | | |
| G46020 | K-FPIP-1300-10BS-5 | 300 | blastocyst handling |
| G26057 | K-FPIP-1600-10BS-5 | 600 | manipulation of oocyte-cumulus complex |
| 10 Vials of 1 (10 Total) | | | |
| G46017 | K-FPIP-1300-10 | 300 | blastocyst handling |

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Customer Service

EU Website: cookmedical.eu
EDI: cookmedical.eu/edi
Distributors: +353 61239240, ssc.distributors@cookmedical.com
Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001633, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Norway: +47 23162968, no.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com

EDI: cookmedical.com/edi.do

Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335

E-mail: orders@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283

E-mail: cau.custserv@cookmedical.com



AI, ESC, IR, OHNS, PI, RH, SUR-A4