


Gain ureteral access with control.

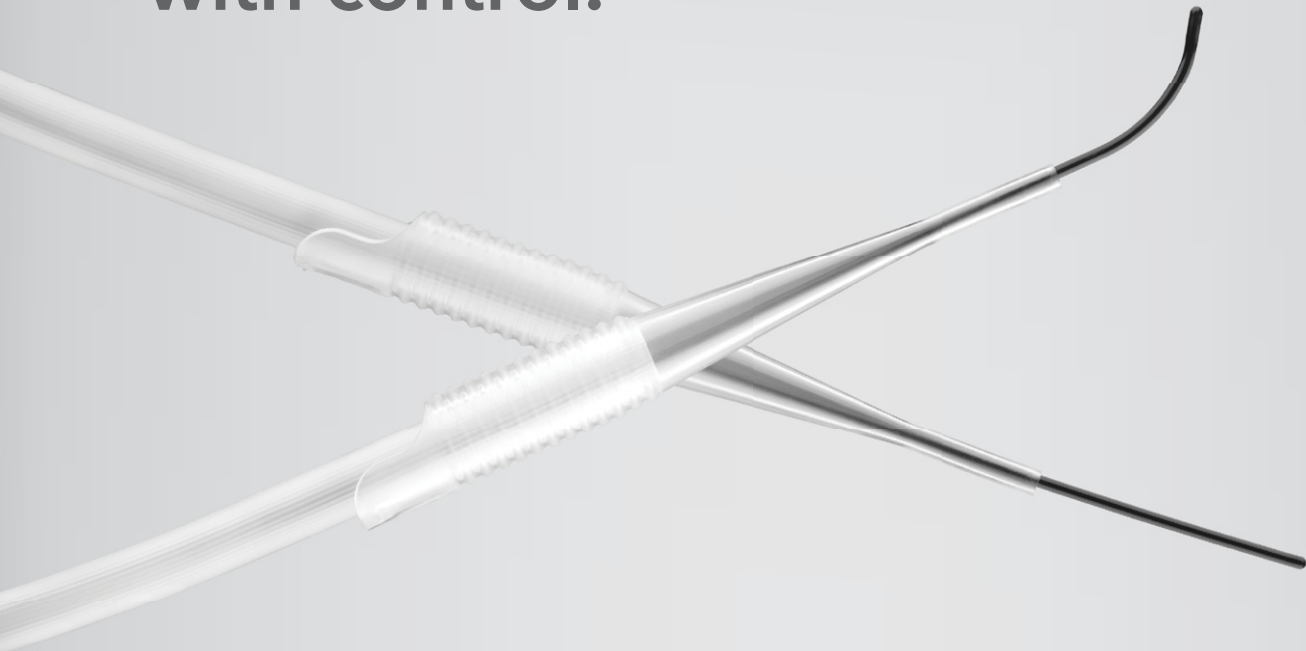

HiWire[®]

NITINOL CORE WIRE GUIDE

Tapered tip

The tapered tip gives you optimized flexibility.


Nitinol core

The nitinol core allows for maximum deflection without kinking.

Hydrophilic coating

The hydrophilic coating reduces friction to provide easier advancement.



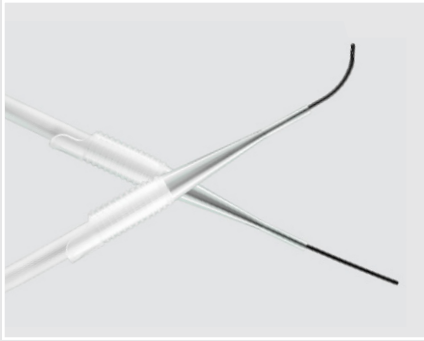

HiWire®

NITINOL CORE WIRE GUIDE

The HiWire® Nitinol Core Wire Guide is used to gain ureteral access in routine procedures or in demanding cases that require increased control. The standard-body and stiff-body configurations allow you to navigate difficult anatomy, straighten tortuous ureters, and place or exchange instruments. The one-to-one torque enables you to maneuver the wire with precision.

Order Number	Reference Part Number	Diameter inch	Length cm	Tip Configuration	Flexible Tip Length cm
Standard-Body Wire Guides					
G30480	HW-025150	.025	150	straight	3
G30476	HW-035150	.035	150	straight	3
G30473	HW-038150	.038	150	straight	3
G30479	HWA-025150	.025	150	angled	3
G30481	HWA-035150	.035	150	angled	3
G30472	HWA-038150	.038	150	angled	3
Stiff-Body Wire Guides					
G30477	HWS-025150	.025	150	straight	3
G30474	HWS-035150	.035	150	straight	3
G30470	HWS-038150	.038	150	straight	3
G30478	HWAS-025150	.025	150	angled	3
G30475	HWAS-035150	.035	150	angled	3
G30471	HWAS-038150	.038	150	angled	3

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Customer Service - Urology

EMEA: EDI - www.cookmedical.com/edi.do
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001633, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Poland: +48 223060159, pl.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com
www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: orders@cookmedical.com

Australia:
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
 E-mail: cau.custserv@cookmedical.com

